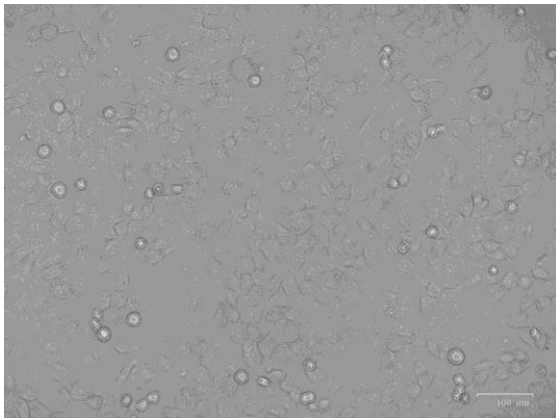


人胃癌细胞 (MKN7 Cells)

细胞货号	T9972	
细胞总数	1 x 10 ⁶ cells	
生长特性	贴壁	
细胞形态	上皮细胞样	
倍增时间	48-72 h	
传代比例	1:2-1:3, 每 2-3 天换液一次	
培养环境	95% 空气, 5% 二氧化碳, 温度 37°C	
培养基	RPMI 1640 (abm Cat. No. TM503) + 10% FBS (abm Cat. No. TM999) + 1% Penicillin/Streptomycin (abm Cat. No. G255), 37.0°C, 5% CO ₂ . 推荐使用: Applied Cell Extracellular Matrix (G422) 包被试剂包被培养皿让细胞更好贴壁。	
冻存条件	70% 完全培养基 + 20%FBS + 10%DMSO	
细胞储存	液氮储存	
细胞培养 注意事项	请参考《细胞培养手册》, 无特殊注意事项。	
STR 数据	- D5S818 : 9,13 - D13S317 : 8,12 - D7S820 : 9,9 - D16S539 : 9,10 - VWA : 17,17 - TH01 : 7,9 - AMEL : X,X - TPOX : 9,9 - CSF1PO : 10,13 - D12S391 : 19,19 - FGA : 23,23 - D2S1338 : 20,23	

电话: 400-833-1885

网址: www.abmGood.com

订购: order.china@china-abm.com

技术: technical.china@abmGood.com



更多细胞详情扫一扫
获取更多信息



下载有问必答app
获取更多信息



	• D21S11 : 29,32.2 • D18S51 : 13,18 • D8S1179 : 15,15 • D3S1358 : 17,17 • D6S1043 : 19,19 • PENTAE : 18,18 • D19S433 : 13,13 • PENTAD : 10,11 • D1S1656 : 11,11
--	---

电话：400-833-1885

网址：www.abmGood.com



更多细胞详情扫一扫
获取更多优惠信息



下载有问必答app
获取更多优惠信息



订购：order.china@china-abm.com

技术：technical.china@abmGood.com